

FORMS, TRACE AND PROCESS

KEEPING TRACK OF THE UNTRACEABLE

Scott Prisbrey, Danny Reyes, Chris Rumbaugh

1



Form Editor <541> Menu

WORKTICKET ENHANCEMENTS

2

WORKTICKET\INVOICE - QR CODE

The Barcode Editor interface includes the following fields:

- Name:
- Length (optional):
- Fixed Text:
- Format:
- Trim Spaces:

➤ Easily add a QR to any Customizable Form

➤ Choose the Barcode option in the ADD TO FORM column

➤ Select QR CODE in the Format dropdown

➤ BENEFIT - QR codes can be scanned from any mobile device using the camera to populate data

➤ Size of QR code is proportional to length of text or amount of data being used

3

WORKTICKET – WKTKA

The WKTKA interface shows a form editor with the following options:

- Form:
- Cursor:
- Height:
- Width:
- Color:
- Group:
- Underline:
- Modify Placement:
- Snap to Grid:

➤ Add customization to the graphic side of the workticket


➤ Requires WKTKA.XMLXF setup from DVI

➤ Enables labs to add data opposite of WKTKB. This can be seen when looking through the tray without removing the workticket

➤ Opportunity to utilize the rotation options

Final Conveyor
USE CAUTION WHEN REMOVING BLOCK!

4

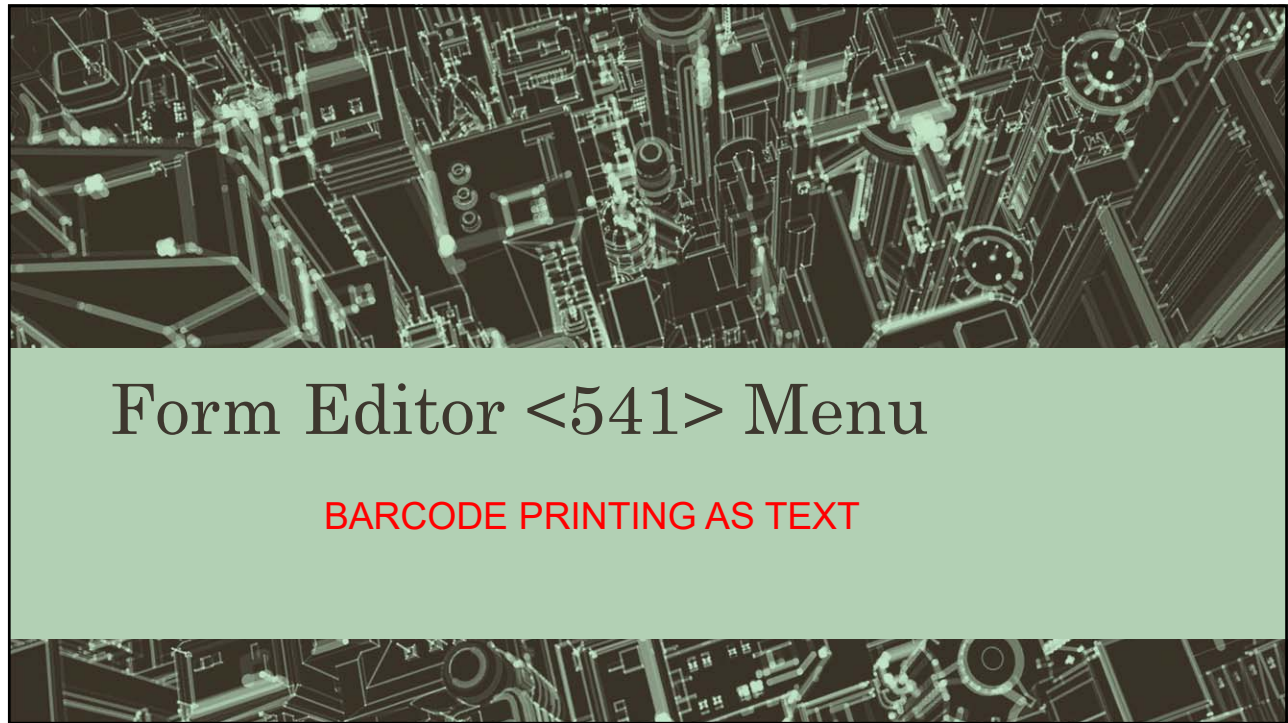

Form Editor <541> Menu

BARCODE PRINTING AS TEXT

5

WORKTICKET – BARCODE PRINTING AS TEXT

- Concatenate data to print on the customizable form
- Select Barcode from the Add To Form column and place it on the template
- Select the data you would like to print and choose to Print As Plain Text from the Format drop down
- Control the font and size using the standard text tools

6

WORKTICKET – OTHER MENTIONABLES

- Ability to Send Data to Front/Back during the form setup
- THJ\$rzsgjGvihw\$Xsgo\$vihw
- Right and Left Lens Size Fields available on Customized Workticket
- New "Frame Received Date" Field on Labels & Worktickets

7



8

WORKTICKET – LABEL PREVIEW

Form

Type InvLbl

Form FTCLBL

Font

Height

Width 0

Color Color

Group

Bold

Italics

Underline

Modify Placement

Right Justify

Snap to Grid

Add to Form

Rectangle

Filled Rectangle

Fixed Text

Barcode

Preview

Cursor

Scr x 401 y 903

Form x 3254.01 y -1542.0

Selection

Scr x y

Form

Label [FTCLBL]

Print

5887 MISSOURI OPTICAL 236308

GD

53/18 04/04/

Enter in-process tray number:

OK Cancel

- Preview the label while editing the setup
- Ability to use in process work to see real-time results

9

Form Editor <541> Menu

DASHBOARD

10

DASHBOARD

REMOTE SETUP <222> MENU

- > Easily designate the Principal remote, which will gather the data
- > Assign Dashboards to other remotes
- > Designated Dashboard will launch when an operator logs into DVI
- > If the Dashboard was active when closedown began it will automatically resume when closedown completes

< > SIGNIFIES A DASHBD CYCLE

Cycles	Dashboard Name	Display Time (secs)
PROMISED	STANDARD	10
SURFBL		5
SURFGE		5
SURFEN		5
PROMISED		10
SURFPO		5
FINBL		5
FINED1		5
PROMISED		10
FINED2		5
FINED3		5
FINED4		5
PROMISED		10
FINF1		5
FINF2		5

11

DASHBOARD

TRACE STATION BY OPERATOR TOTAL


- The ability to track daily total numbers by operator

Operator	Count
BB	0
AH	3
JB	8

12

DASHBOARD

TRACE STATION BY OPERATOR HOURLY SHIFT

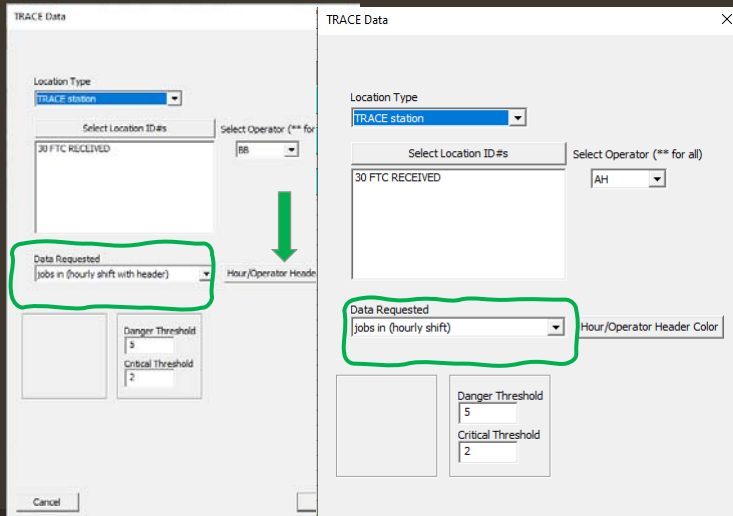

- Track trace station hourly totals for a 9-hour shift by operator
- Choose the 1st operator and select Hourly Shift With Headers
- Control the Hour/Operator Header colors separately from the threshold
- Insert other operators without the headers

13

DASHBOARD


TRACE STATION BY OPERATOR HOURLY SHIFT

- THE ABILITY TO TRACK HOURLY NUMBERS BY OPERATOR OVER A 9 HOUR SHIFT

FTC RECEIVED

	8:00	9:00	10:00	11:00	12:00	13:00	14:00	15:00	16:00
BB	0	0	0	0	0	0	0	0	12
AH	0	0	0	0	0	0	3	0	3
JB	0	0	0	0	0	0	8	0	0

14

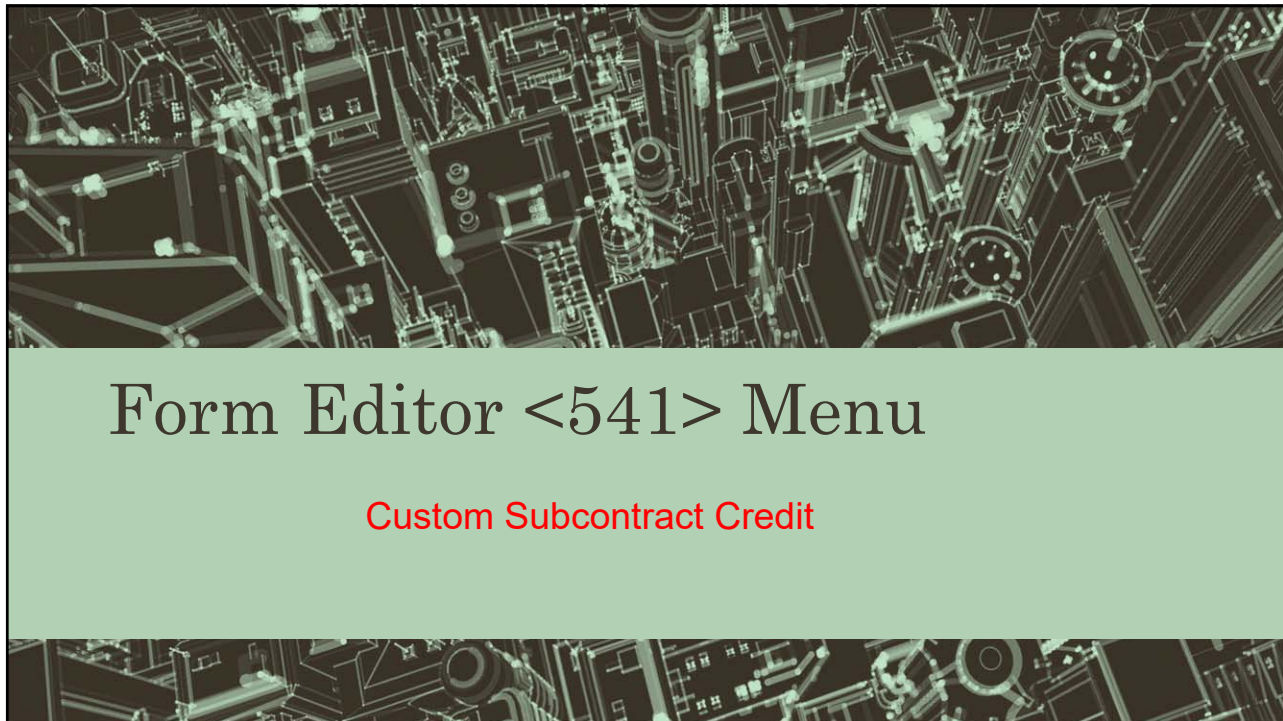

15

Custom Subcontract Credit

- Subcontract credits now customizable in <541>

Form

Type: Subcontract Credit

Form: CSUBFORM

CSUBCUST

CSUBFORM

Height: [dropdown]

16



17

Station Analyzer


- Ability to look at tray activity by machine group and specific stations within a group.
- Ability to show job details such as styles, materials, edge types, coatings and frames.
- Ability to show elapsed times and averages at a station or between stations by group and individual trays.
- Exporting of report group and detail data is available.

18

Station Analyzer

Screen Layout


- 1 - Scan Configuration
Set parameters for job inclusion, station type or location, job group detail, cycle time calculation, and report time span .
- 2 - Job/Tray Group Output
Report output subtotals, including usage counts and timing averages. Defined by group selection
- 3 - Job/Tray Detail Output
Job detail by group

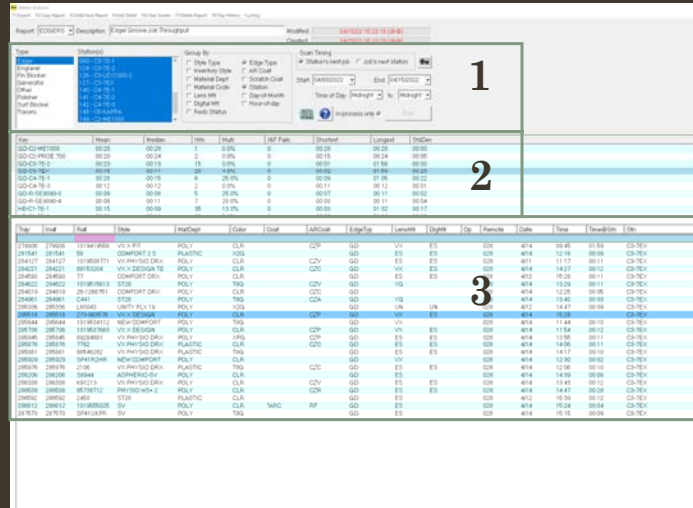

19

Station Analyzer

- Stations types (Edger, Generator, etc.) are defined in a machines function setup in <487>. Stations without a machinery function reside in the “other” category.
- Numerous options for grouping jobs

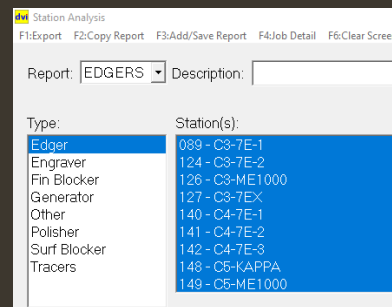

20

Station Analyzer

- Report on tray activity by days, hours or shifts and whether to report times by "Station's next job" (JPH) or "Job's next station"

Scan Timing

Station's next job Job's next station

Start: 04/08/2022 End: 04/15/2022

Time of Day: Midnight to Midnight

In-process only Run

21

Station Analyzer

- Selecting a group in the group panel will show job level data in the detail panel

Key	Mean	Median	Stk	Multi	INF Fails	Shortest	Longest	StdDev
GD-C2-ME1000	00:28	00:28	1	0.0%	0	00:28	00:28	00:00
GD-C2-PRDE 700	00:20	00:24	2	0.0%	0	00:15	00:24	00:05
GD-C3-7EX	00:15	00:11	23	4.8%	0	00:02	01:59	00:28
GD-C4-7C-3	00:12	00:12	2	0.0%	0	00:11	00:11	00:01
GD-R-SE9090-3	00:09	00:08	5	25.0%	0	00:07	00:02	00:02
GD-R-SE9090-4	00:08	00:11	7	20.0%	0	00:03	00:04	00:04
HF-C1-7E1	00:16	00:09	35	13.3%	0	00:03	00:17	00:17

Tray	Intr	Out	Style	MatDept	Color	Coat	ARCcoat	EdgeTyp	LensMtr	DgMtr	Op	Remote	Date	Time	Time@Stn	Stn
279028	279036	1019419558	VX-X FIT	POLY	CLR		C2P	GD	VX	ES	028	414	09:46	01:59		C3-TEX
281541	281541	58	COMFORT 2 S	PLASTIC	X2G			GD	ES	ES	028	414	12:16	00:09		C3-TEX
284127	284127	1019509771	VX PHYSIO DRX	POLY	CLR		C2V	GD	ES	ES	028	411	11:17	00:11		C3-TEX
284221	284221	80165004	VX DESIGN TE	POLY	CLR		C2C	GD	VX	ES	028	414	14:27	00:12		C3-TEX
284593	284593	77	COMFORT DRX	POLY	CLR			GD	ES	ES	030	412	15:23	00:11		C3-TEX
284622	284622	1019515613	ST2	POLY	TRG		C2V	GD	YG		028	414	13:29	00:11		C3-TEX
284819	284819	20-126761	COMFORT DRX	POLY	CLR		C2C	GD	VX		028	414	12:25	00:05		C3-TEX
284961	284961	0441	ST28	POLY	TRG		G2A	GD	YG	UN	028	414	13:40	00:05		C3-TEX
285336	285336	L50043	UNITY PLX 19	POLY	X2G			GD	UN	UN	028	412	14:47	00:09		C3-TEX
285818	285818	27906826	VX-X DESIGN	POLY	CLR		C2P	GD	VX	ES	028	414	15:28			C3-TEX
285844	285844	1019534112	NEW COMFORT	POLY	TRG			GD	VX	ES	028	414	11:44	00:10		C3-TEX
285796	285796	1019537963	VX-X DESIGN	POLY	CLR		C2P	GD	VX	ES	028	414	11:54	00:12		C3-TEX
285846	285846	80264981	VX PHYSIO DRX	POLY	XPS			GD	ES	ES	028	414	13:55	00:11		C3-TEX
285876	285876	7782	VX PHYSIO DRX	PLASTIC	CLR		C2C	GD	ES	ES	028	414	14:06	00:11		C3-TEX
285881	285881	88584232	VX PHYSIO DRX	PLASTIC	TRG			GD	ES	ES	028	414	14:17	00:10		C3-TEX
285929	285929	SP41762R	NEW COMFORT	POLY	CLR			GD	VX		028	414	12:30	00:02		C3-TEX
285976	285976	2196	VX PHYSIO DRX	PLASTIC	TRG		C2C	GD	ES	ES	028	414	12:06	00:10		C3-TEX
286296	286296	33944	ASPHERIC SV	POLY	CLR			GD	ES	ES	028	414	14:59	00:08		C3-TEX
286338	286338	K91213	VX PHYSIO DRX	POLY	CLR		C2V	GD	ES	ES	028	414	13:43	00:12		C3-TEX
286838	286838	86798712	PHYSIO W5+ 2	POLY	CLR		C2B	GD	ES	ES	028	414	14:47	00:28		C3-TEX
286952	286952	2450	ST28	FLASTIC	CLR			GD	ES		028	412	15:59	00:12		C3-TEX
286612	286612	1019555025	SV	POLY	CLR	*ARC	RF	GD	ES		028	414	15:24	00:04		C3-TEX
287673	287673	SP41UKFR	SV	POLY	TRG			GD	ES		028	414	15:16	00:09		C3-TEX

22

Station Analyzer

Options

- Rt-Click on Job Detail header to select display columns
- F4 Job Detail (Rx job screen)
- F8 Tray history

Job	Mean	Median
GD-CA-7E-1	00:29	00:15
GD-CA-7E-9	00:12	00:12
GD-R-SE0002	00:06	00:09
GD-R-SE0003	00:09	00:08
GD-R-SE0004	00:08	00:11
HB-CI-7E-1	00:15	00:09
HB-CI-7E-2	00:18	00:11
HB-CI-7E-3	00:14	00:08
HB-CI-7E-4	00:18	00:09

Tray	Tray#	Ref#	Style	Color	Coeff	AP
295569	295569	00-1214006	SV			
295526	295526	0026	SV			
279214	279214	101948709	SV			
280174	280174	1019422416	COMFORT DR			
281186	281186	1019482482	SV PLYD DR			
284803	284803	1019821092	QUALTRON			
285021	285021	1019535660	AUTO 10 1 9			
285906	285906	1019543737	ST28			
286029	286029	1019555015	SV			
287098	287098	1019560098	SV			
287676	287676	1019591863	SV			
288033	288033	1019596051	SV			
282560	282560	1130	ST28			
282726	282726	117108	SV			
285440	285440	157201	SV			
6216	160586	2288177	SV			
272761	272761	27213	SV			
272761	272761	27213	SV			
272761	272761	27213	SV			
283940	283940	9576	SV			
281962	281962	37967	SV			
284126	284126	38-1288181	IGDAL			
281951	281951	40-1288318	COMFORT DR			
284299	284299	46-1288354	COMFORT DR			
281967	281967	461606	SV			
283024	283024	41-1287281	SV			
283024	283024	41-1287281	SV			
282922	282922	52-1288764	COMFORT DR			
286791	286791	52-1288616	SV			
283077	283077	5384	SV			

Station Analysis

F1:Export F2:Copy Report F3:Add/Save Report **F4:Job Detail** F6:Clear Screen F7>Delete Report **F8:Tray History**

Report: TRACERS Description:

Type: Edger Engraver Fin Blocker Generator

Station(s): 072 - R-SE000-1 073 - R-SE000-2 074 - R-SE000-3 075 - R-SE000-4

Group By: Style Type Inventory Style Material Dept

23

Station Analyzer

- Save report setups

Station Analysis

F1:Export F2:Copy Report **F3:Add/Save Report** F4:Job Detail F6:Clear Screen F7>Delete

Report: TRACERS Description:


Type: **TRACERS** Edger Engraver Fin Blocker Generator Other Polisher Surf Blocker Tracers

Station(s): 001 - TRACER #1 002 - TRACER #2 003 - TRACER #3 004 - TRACER #4 005 - TRACER #5 138 - TRACER #6

24

Station Analyzer


- Cycling can run a saved report automatically

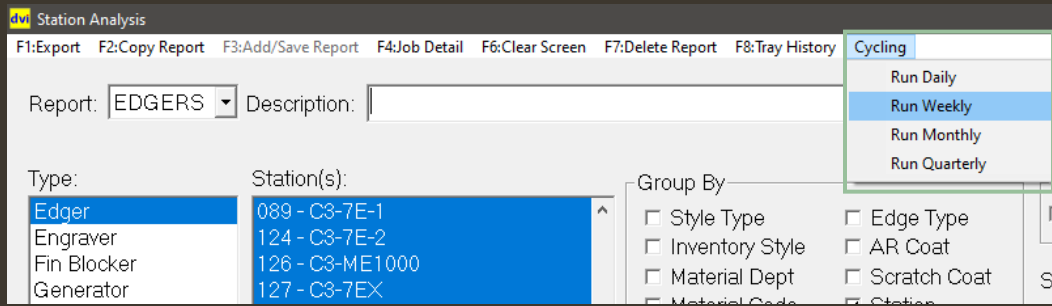

25

Station Analyzer

- F1 to export report (*.csv) to view and work data in Excel

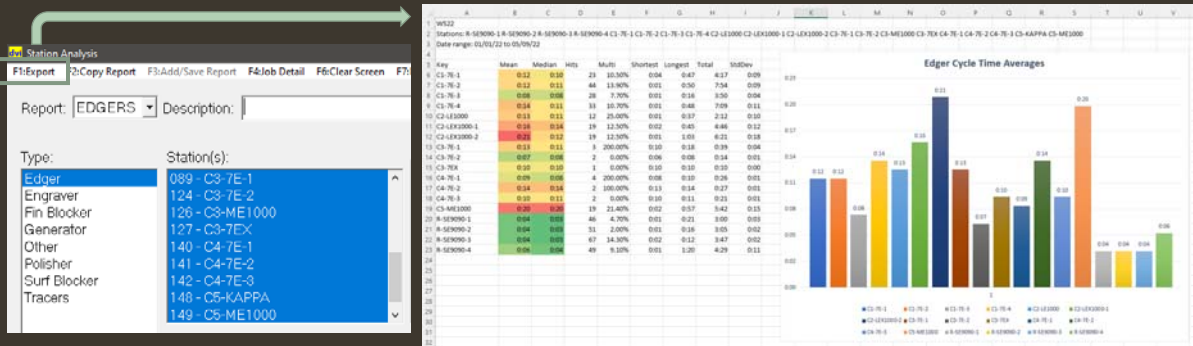

26



New Label Printing Interface

Printing Labels

27

Print By Account

- <2365> in the menus
- Dropdown allows you to select the desired form
- Running list of entries
- Set number of copies

Digital Vision, Inc 05/03/2022
 6805 SE Milwaukie Ave
 Portland, OR 97202


YVETTE CANADA OPTICAL
 6805 SE MILWAUKEE
 PORTLAND CN 02888

6

28

Routing Features


- <2259> in the menus
- Dropdown with filter logic display
- Running list of entries
- Log of scans
- Scannable by tray or invoice

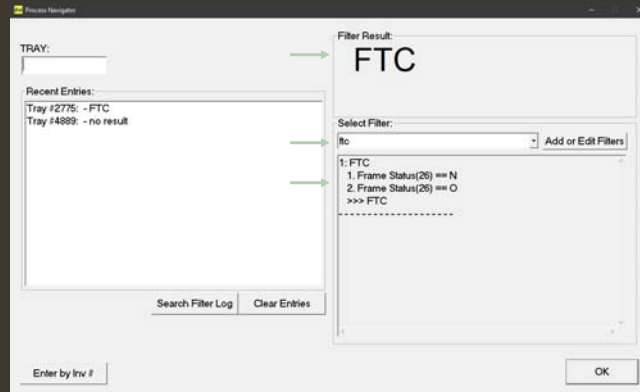

29

Search Prior Scans

- Click "Search Filter Log"
- Enter an invoice number

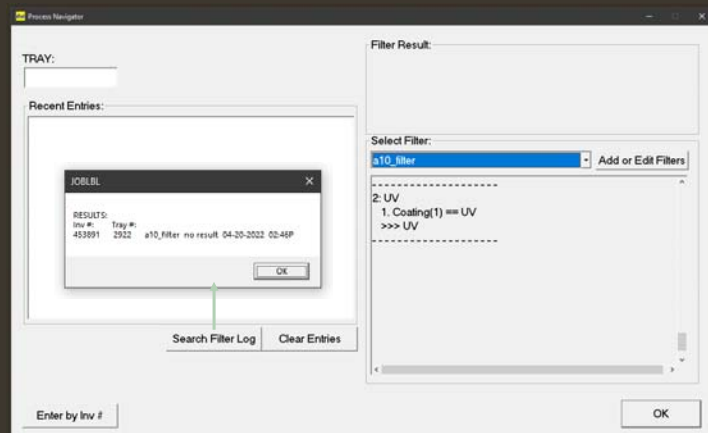

30



31